# MASS CALIBRATION AND RELATIVE HUMIDITY COMPENSATION REQUIREMENTS FOR OPTICAL PORTABLE PARTICULATE MATTER MONITORS: THE IMPASHS (IMPACT OF SMOKE-FREE POLICIES IN EU MEMBER STATES) WP2 PRELIMINARY RESULTS

Ruprecht AA<sup>1,2</sup>, De Marco C<sup>1</sup>, Boffi R<sup>1</sup>, Mazza R<sup>1</sup>, Lopez MJ<sup>3</sup>, Moshammer H<sup>4</sup>,

Dautzenberg B<sup>5</sup>, Clancy L<sup>6</sup>, Precioso J<sup>7</sup>, & Invernizzi G<sup>1,2</sup>

<sup>1</sup> Fondazione IRCCS Istituto Nazionale Tumori di Milan Italy, <sup>2</sup> SIMG-Società Italiana di Medicina Generale Italy, <sup>3</sup> Agencia de Salut Publica de Barcelona Spain, <sup>4</sup> Medical University of Vienna Austria, <sup>5</sup> OFT Paris France, <sup>6</sup> Research Institute for a Tobacco Free Society Dublin, Ireland, <sup>7</sup> Instituto de Educação e Psicologia. Universidade do Minho, Braga, Portugal.

info@tecanalysis.it

#### **Abstract**

#### Introduction

Better knowledge of particulate matter (PM) concentrations needs portable, reliable, user friendly, low cost, real time mass analyzers of PM2.5 and PM10.

Optical Particle Counters (OPC) measuring mass have manufacturer calibration specific gravity "K" factor referred to polystyrene latex particles which are completely different than those of the real world, therefore they require specific calibrations. Measurements are also subject to Relative Humidity (RH) heavy interference.

### **Objectives**

To evaluate, within the IMPASHS WP2 Project, the performance of four different OPC's in Environmental Tobacco Smoke and background urban pollution and to find the new "K" factors using one Model BAM-1020 with certificate n. EPQM-0798-122 for comparison.

#### Methodology

All instruments have been operating in parallel measuring PM2.5 generated by cigarettes (ETS) indoor and by urban pollution outdoor and the data were replicated three times.

#### Results

"K" factors were widely different between manufacturer's model, instrument serial numbers, ETS and urban pollution, ranging from 0.5 to 2.27. Correlation with BAM-1020 was ranging from 0.7500 to 0.9800 and Student t test from 0.3000 to 0.9500. Relative Humidity interference resulted mathematically compensable up to 75 % RH, but above becomes uncontrollable and sample drying becomes compulsory.

#### Conclusion

OPC's are very reliable and accurate but need specific calibration and special care in handling and elaboration of the measurements.

### Introduction


Measurement of indoor/outdoor particulate matter (PM) pollution in real time can be satisfactorily and reliably performed using Optical Particle Counters (OPC), if properly calibrated according to the procedure described in the next paragraphs. The OPC's principle of operation is the nephelometric measurement based on the light scattering of airborne particles. The sample of air is normally drawn into the light scattering sensor (nephelometer) with a flow controlled rotary vane pump. The nephelometer is constituted by a light source (low power laser diode), scattered light collection optics and a photo detector circuit. The flow path for the air crosses the path of the laser diode. When the air is clean (absence of airborne particulate), the laser diode light is extinguished in a light trap, but when the air is containing airborne particulate, there is scattering of light which is collected and measured.

The output signals from the detector are empirically proportional to the number and size of particles and are elaborated to present the data expressed in number of particles per liter of selected classes of diameters.

However the certified particulate matter measurement systems are based on gravimetric methods and the concentrations are expressed in micrograms per cubic meter ( $\mu g/m^3$ ) therefore it is necessary to convert the number of particles per liter to mass. Unfortunately all light scattering devices have inherent difficulties when converting light scatters to mass. Index of refraction and mean particle diameter can affect the amount of light scattered from the same amount of mass<sup>6</sup>. A mathematical equation is applied to the number/size of particles measured by the OPC to obtain first the particle total volume on which a "K" factor is applied to adjust for the specific gravity of the PM to be measured. The simplest solution is to compare the OPC's mass measurements made with default factory "K"=1.000 for a set period of time with the mass measured by a gravimetric system over the same period of time. Comparing the concentrations will yield a "K" that can be applied to all OPC measurements performed by the same calibrated analyzer.

Another factor which greatly influences the mass measurements of the OPC's is the interference of the relative humidity (RH). Sioutas *et al.* in 2000, and Chakrabarti *et al.* in 2004 showed that RH could drastically affect the mass concentrations measurements of one OPC model as demonstrated in Fig. 1 in which the desired concentration ratio is 1 and, at RH values greater than 50 %, this ratio begins to increase due to particle aggregation (particle size increase as water is absorbed).

# Humidity impacts on pDR

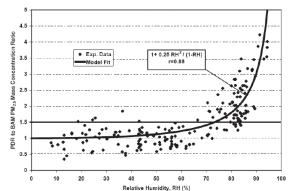

Thermo pDR PM 2.5 mass readings compared to a reference mass monitor. Readings become non linear at RH above  $\sim$ 50-60%.

B. Chakrabarti et al., 2004

Fig. 1

To prevent this inaccuracy it is necessary to provide the OPC measurements with a simultaneous RH measurements and compensate mathematically. However above 80/85 % RH, the interference become so high that inaccuracy in RH measurements brings to unacceptable errors in the concentrations. To obviate to this errors and to permit the use of OPC's also when RH is greater than 80 % it become necessary to reduce the RH to the light scatter chamber to about 50 % either by heating the sample or by drying it by means of *permapure* dryer<sup>TM</sup>.

## **Objectives**

To evaluate, within the IMPASHS WP2 Project, the performance of four different OPC's in the measurements of Environmental Tobacco Smoke (ETS) and background urban pollution and to find the new "K" factors for PM<sub>2.5</sub> using as automatic mass reference the measurements of one Model BAM-1020 with certificate of equivalence

## Methodology

The method to measure PM environmental pollution used in monitoring stations all over the world is the gravimetric (Federal Reference Method FRM or the European equivalent) and the data are generally presented as mass of PM<sub>2.5</sub> and PM<sub>10</sub>. The system consists of equipment designed to accumulate the PM on a pre-weighted filter through which a known air flow is sampled for a programmable time. At the end of the sampling time the filter is removed and accurately dried and weighted. The difference in weight is the amount of PM deposited on the filter which, divided for the number of cubic meters of air passed through the filter, gives the mass expressed in micrograms per cubic meter (μg/m<sup>3</sup>). However this method requires skilled operators and expensive laboratory equipment, is time consuming and deliver the results of the measurements only with a delay of some days. To overcome this inconvenient, automatic sampling and measurement methods have been developed which can deliver the concentrations every 1 or 2 hour's. These methods are the Beta Attenuation Monitors (BAM) and Tampered Element Oscillating Balance (TEOM). Several manufacturers have developed analyzers based on these principles of operation and many have been designed as equivalent to the gravimetric by the U.S.A. Environmental Protection Agency. These analyzers are commonly used in most of the monitoring stations and deliver the concentrations in hourly averages informing in almost real time the pollution level. The PM<sub>2.5</sub> measurements of the OPC models described in Appendix A are compared with the same PM<sub>2.5</sub> measured by one Beta Attenuation Monitor model BAM-1020, equipped with standard PM<sub>2.5</sub> inlet. The BAM-1020, is manufactured by Metone Instruments Inc. and is designated as equivalent method by U.S. E.P.A. n.° EQPM-0798-122 and certified by the German T.Ü.V. (: 936/21205333/A Köln, 06.12.2006).

**Site informations:** the generation of ETS aerosol and the measurements have been performed in the Tobacco Control Unit Research Laboratory of the Fondazione IRCCS Istituto Nazionale dei Tumori, located in Milan, in a room of about 45 m<sup>3</sup> without air conditioning and with about 0.3/0.4 air exchange per hours (ach). To assure the maximum mixing factor (PM concentrations of the same value in every point of the room), one fan of about 1,500 m<sup>3</sup>/hour was always in operation during the ETS generation and measurements.

### **Results**

## **RH** Interference compensation

During all tests the RH was measured and all raw data have been mathematically compensated according the Sioutas tests using the following equation:


$$C_{comp} = C_{meas} / ((1 + (RH/100)7) \times 3.72)$$

Where:  $C_{comp}$  = concentration compensated for RH interference

 $C_{meas}$  = concentration measured

**RH** = Relative Humidity measured

Fig.2 shows the interference compensation curve and Fig. 3 the linear regression analysis.

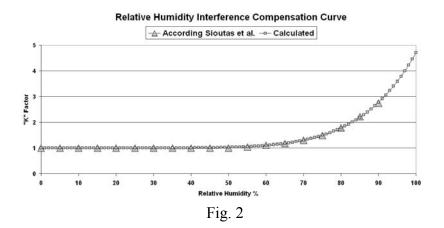

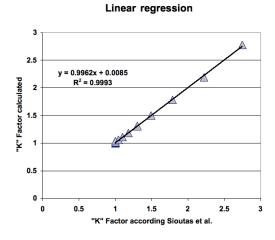


Fig. 3

The accuracy of the RH interference correction using the above equation has been experimentally tested measuring continuously the outdoor PM<sub>2.5</sub> concentrations over a wintertime period for 7 days using two OPC's: the model e-sampler and the model Aerocet 531 (same instruments used for the gravimetric calibrations) operating in parallel.

The e-sampler is normally equipped with a sampler heater control system driven by a RH sensor installed on the inlet of the laser chamber and the RH set point is programmable. In this case the RH set point was set at 40 %. The Aerocet is equipped only with a Temperature and RH measurement.

The results are represented in Fig. 4: the RH changed from a minimum of about 30 to a maximum of 98 % with no interference on the e-sampler indications (red line) but heavily interfering with the Aerocet 531 indication by increasing the value of a factor of 3 with RH > 80/85%. The thin black line is showing the uncorrected Aerocet 531 measurements. When the equation was applied, in excel, to the Aerocet 531 raw data, the correction was very effective. Correlation analysis and Student t test between e-sampler and Aerocet 531 indicate the accuracy of the correction giving the following results:

Without correction: correlation = 0.8511 and p = 7.7E-49

After correction: correlation = 0.9765 and p = 0.6859


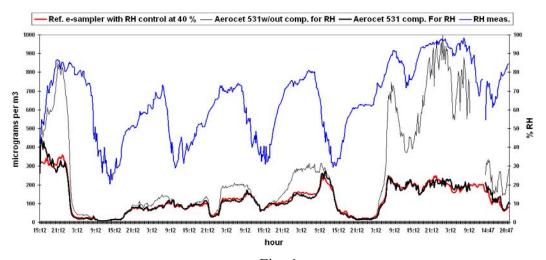


Fig. 4


## "K" factor corrections

All Analyzer measurements were made using the default "K" factor of 1.000 and the comparison yields to the results listed below for each instrument. Comparison is made on a hourly basis since the BAM-1020 supply the data as the average. This operation was made in Excel calculating the hourly average of the Analyzers with their internal clock synchronized with the BAM-1020 clock before starting the sampling. For each set of data, the correlation and the "p" of Student t test is calculated and reported to evaluate the accuracy and precision of the new "K" factor resulting from the calibration.

An example of the graphs and the regression analysis is shown in Fig. 5 and 6:

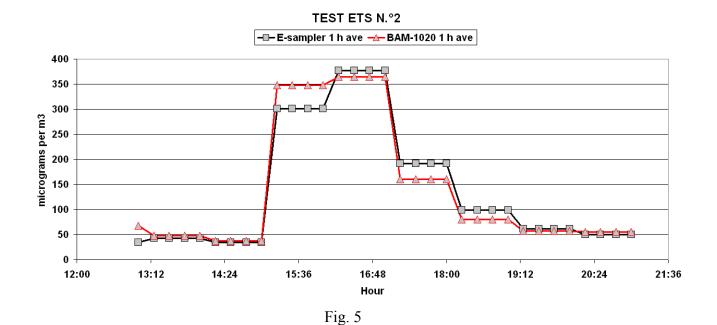

## ESAMPLER/BAM-1020 TEST N.°2 REGRESSION ANALYSIS

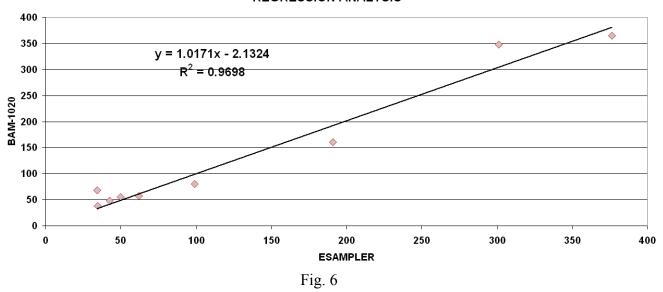

The following tables show the detailed results for each test and each OPC.

# 1.- Model Aerocet 531, Metone Instruments Inc. serial # E-1871

## Tobacco Control Unit, National Cancer Institute, Milan, Italy

III° test	BAM	Aerocet	Correlation	P test	K factor
ETS	128.0	129.2	0.9147	0.2361	11.5
URBAN	49.0	49.0	0.7943	0.9769	8.1
II° test	BAM	Aerocet	Correlation	P test	K factor
ETS	179.8	179.8	0.9940	0.9767	10.2
URBAN	57.3	57.4	0.8574	0.7661	14.3
I° test	BAM	Aerocet	Correlation	P test	K factor
ETS	201.5	201.6	0.9912	0.9478	12.2
URBAN	52.0	52.0	0.9032	0.9871	16.9
Average	Correlation	P test	K factor (SD)		
ETS	0.9666	0.7202	11.3 (1.0)	1	
URBAN	0.8516	0.9100	13.1 (4.5)	1	

# 2.- Model Aerocet 531, Metone Instruments Inc. serial # 8551 $\,$

# Research Institute for a Tobacco Free Society, Dublin, Ireland

BAM	Aerocet	Correlation	P test	K factor
189.1	186.8	0.9821	0.9310	10.2
18.9	19.0	0.9665	0.5470	6.6
BAM	Aerocet	Correlation	P test	K factor
133.2	133.4	0.9915	0.7079	8.1
na	na	na	na	na
BAM	Aerocet	Correlation	P test	K factor
86.3	86.6	0.9871	0.8629	9.2
na	na	na	na	na
Correlation	P test	K factor (SD)		
0.9869	0.8340	8.1 (1.05)		
na	na	6.6 (na)		
	189.1  18.9 <b>BAM</b> 133.2  na <b>BAM</b> 86.3  na <b>Correlation</b> 0.9869	189.1       186.8         18.9       19.0         BAM       Aerocet         133.2       133.4         na       na         BAM       Aerocet         86.3       86.6         na       na         Correlation       P test         0.9869       0.8340	189.1       186.8       0.9821         18.9       19.0       0.9665         BAM       Aerocet       Correlation         133.2       133.4       0.9915         na       na       na         BAM       Aerocet       Correlation         86.3       86.6       0.9871         na       na       na         Correlation       P test       K factor (SD)         0.9869       0.8340       8.1 (1.05)	189.1       186.8       0.9821       0.9310         18.9       19.0       0.9665       0.5470         BAM       Aerocet       Correlation       P test         133.2       133.4       0.9915       0.7079         na       na       na       na         BAM       Aerocet       Correlation       P test         86.3       86.6       0.9871       0.8629         na       na       na       na         Correlation       P test       K factor (SD)         0.9869       0.8340       8.1 (1.05)

# 3.- Model Aerocet 531, Metone Instruments Inc. serial # F-8557

# **OFT, Paris, France**

				K factor
189.1	189.1	0.9818	0.1924	8.3
18.9	18.9	0.9777	0.8638	5.5
BAM	Aerocet	Correlation	P test	K factor
109.3	109.4	0.9904	0.9278	6.3
na	na	na	na	na
BAM	Aerocet	Correlation	P test	K factor
86.3	86.7	0.9974	0.5606	7.0
na	na	na	na	na
Correlation	P test	K factor (SD)		
0.9899	0.5603	8.1 (1.01)		
na	na	5.5 (na)		
	18.9  BAM  109.3  na  BAM  86.3  na  Correlation  0.9899	18.9       18.9         BAM       Aerocet         109.3       109.4         na       na         BAM       Aerocet         86.3       86.7         na       na         Correlation       P test         0.9899       0.5603	18.9   18.9   0.9777	18.9   18.9   0.9777   0.8638

# 4.-Model 1.180 Portable Aerosol Spectrometer, Grimm serial #8F070042

# Inst. Umwelthygiene, ZPH, Wien, Austria

I° test	BAM	Grimm	Correlation	P test	K factor
ETS	253.3	253.0	0.9990	0.9625	1.26
URBAN	51.7	51.5	0.8423	0.6725	1.30
II° test	BAM	Grimm	Correlation	P test	K factor
ETS	198.2	199.2	0.9781	0.9468	1.22
URBAN	56.9	54.4	0.9733	0.4936	1.30
III° test	BAM	Grimm	Correlation	P test	K factor
ETS	151.6	152.7	0.9969	0.7579	1.10
URBAN	49.6	49.6	0.9391	0.9308	1.26
Average	Correlation	p test	K Factor (SD)		
ETS	0.9914	0.8890	1.2 (0.08)	-	
URBAN	0.9182	0.6990	1.3 (0.02)	-	

# 5.- Model SidePak TSI Serial# 10805037

## Servei d'Avaluaciò i Metodes d'Intervenciò

# Agencia de Salut Publica de Barcelona, Barcelona, Spain

III° test	BAM	S. 10805037	Correlation	P test	K factor
ETS	135.8	134.7	0.9980	0.7397	0.45
URBAN	49.1	49.1	0.9304	0.9763	0.51
II° test	BAM	S. 10805037	Correlation	P test	K factor
ETS	198.2	49.1	0.9304	0.9763	0.51
URBAN	58.0	57.6	0.9756	0.6299	0.47
I° test	BAM	S. 10805037	Correlation	P test	K factor
ETS	196.5	196.3	0.9947	0.9880	0.56
URBAN	52.8	52.8	0.8990	0.9700	0.52
Average	Correlation	P test	K factor (SD)		
ETS	0.9744	0.9013	0.51 (0.05)		
URBAN	0.9350	0.8587	0.50 (0.02)		

# 6.- Model SidePak TSI Serial# 10805044

# Servei d'Avaluaciò i Metodes d'Intervenciò

# Agencia de Salut Publica de Barcelona, Barcelona, Spain

III° test	BAM	S. 10805044	Correlation	P test	K factor
ETS	135.8	135.9	0.9981	0.9727	0.45
URBAN	49.1	49.5	0.9263	0.6171	0.48
II° test	BAM	S. 10805044	Correlation	P test	K factor
ETS	157.6	157.2	0.9884	0.9677	0.52
URBAN	47.2	47.6	0.3453	0.6901	0.48
I° test	BAM	S. 10805044	Correlation	P test	K factor
ETS	176.0	174.9	0.9955	0.9020	0.55
URBAN	52.7	52.4	0.9184	0.8267	0.55
Average	Correlation	P test	K factor (SD)		
ETS	0.9940	0.9475	0.51 (0.05)	-	
URBAN	0.7300	0.7113	0.50 (0.04)		

## 7.- Model e-sampler, Metone Instruments Inc. serial # G3427

## Tobacco Control Unit, National Cancer Institute, Milan, Italy

III° test	BAM	e-sampler	Correlation	P test	K factor
ETS	224.0	224.6	0.8621	0.4453	1.77
URBAN	52.4	52.0	0.9659	0.4390	1.85
II° test	BAM	e-sampler	Correlation	P test	K factor
ETS	141.6	141.3	0.9884	0.9070	1.77
URBAN	56.9	56.6	0.9884	0.9070	1.66
I° test	BAM	e-sampler	Correlation	P test	K factor
ETS	123.3	124.2	0.9568	0.5634	1.51
URBAN	49.2	49.1	0.9568	0.8593	1.24
Average	Correlation	P test	K factor (SD)		
ETS	0.9358	0.6385	1.68 (0.1)		
URBAN	0.9704	0.7351	1.58 (0.3)		

## **Summary of "K" factors**

	Aerocet	Aerocet	Aerocet	Grimm	SidePak	SidePak	e-sampler
	# 1	# 2	# 3	# 4	# 5	# 6	# <b>7</b>
ETS	11.3 (1.0)	8.1 (1.05)	8.1 (1.01)	1.2 (0.08)	0.51 (0.05)	0.51 (0.05)	1.68 (0.1)
URBAN	13.1 (4.5)	6.6 na	5.5 na	1.3 (0.02)	0.50 (0.02)	0.50 (0.04)	1.58 (0.3)

## **Comments**

"K" factors were found to vary noticeably not only between manufacturer's model but also between instruments of the same manufacturer with different serial numbers.

Differences in the "K" factors between ETS and urban pollution have also been found and must be taken in consideration during the measurements campaigns.

Long term measurements in urban pollution (7 days) showed changes in the optimal "K" factor suggesting that there may be changes in urban pollution composition affecting it sensibly, but the evaluation of the amount of these changes needs more investigations.

These facts are confirming the very well known phenomena of optical characteristics modifications caused by changes in the physical/chemical composition and morphology of the PM.

These factors together with the modification in specific gravity are consequently changing the final mass measurement whenever the typology of PM is changing.

In the test conditions, the "K" factors were ranging from 0.5 to 13.1. But despite this wide range, the accuracy and precision for a given typology of PM measured with all OPC's compared with the BAM-1020 measurements were good, ranging from about 0.3500 to 0.9800 and Student t test from about 0.3000 to 0.9500. Also the repeatability of the measurements of all analyzers when used with the same PM characteristics is acceptable.

#### Conclusion

All OPC's can be considered suitable for the measurements of ETS PM provided the above "K" factors are applied to the raw data. Regarding urban pollution it must be advised that there may be differences in case the physical/chemical composition and morphology of the PM of the site where the measurements are made differs from those of Milan, where the calibration have been made. In this case it is suggested to proceed locally with a further calibration to determine the new factor.

Relative Humidity interference can be compensated mathematically and with acceptable accuracy up to 80/85 % RH, but above these values the error introduced by the RH measurement (normally  $\pm 5$  % RH) becomes unacceptable and sample drying becomes compulsory.

OPC's are extremely simple and user friendly to use, maybe too simple: just push a button and they go and in a few minutes they produce the data. They are also very reliable and accurate. But the correct interpretation of the data is not so simple; manufacturer factory calibration data cannot be applied and they absolutely need

specific calibration and special care in handling and elaboration of the measurements applying the procedures above described.

### References

1.- Repace JL, Lowrey AH. Indoor air pollution, tobacco smoke and public health. *Science* 

1980;208:464–72.

- 2.- Invernizzi G., Ruprecht A.A., De Marco C., Mazza R., Rossetti E., Boffi R.. Real-time measurement of indoor particulate matter originating from environmental tobacco smoke: a pilot study. *Epidemiol Prev* 2002;26:30–4.
- 3.- Boffi R., Ruprecht A.A., Mazza R., Ketzel M., Invernizzi G.. A day at the European Respiratory Society Congress: passive smoking influences both outdoor and indoor air quality. *Eur Resp. J*, 2005, 4: 862-863.
- 4.- Repace J. Measurements of outdoor air pollution from secondhand smoke on the UMBC campus. http://www.repace.com/pdf/outdoorair.pdf. Date last accessed: December 2005. Date last updated: November 2005.
- 5.- Invernizzi G. Ruprecht A.A., De Marco C., Mazza R., Lopez M. J., Gorini G., Boffi R... Outdoor Environmental Tobacco Smoke: its contribution to outdoor ambient pollution in a selected area of the historical city center of Milan. *In publication*
- 6.- Jenkins R.A., Ilgner R.H., Tomkins B.A., Peters D.W. Development and Application of Protocols for the Determination of Response of Real-Time Particle Monitors to Common Indoor Aerosols *Air & Waste Manage*. *Assoc.* 2004 54:229–241
- 7.- Chung, A.; Chang, D.P.Y.; Kleeman, M.J.; Perry, K.D.; Cahill, T.A.; Dutcher, D....; Stroud, K. Comparison of Real-Time Instruments Used to Monitor Airborne Particulate Matter; *J. Air & Waste Manage*. *Assoc.* 2001, *51*, 109-120.
- 8.- Ott, W.; Switzer, P.; Robinson, J. Particle Concentrations inside a Tavern before and after Prohibition of Smoking—Evaluating the Performance of an Indoor Air-Quality Model; *J. Air & Waste Manage. Assoc.* 1996, *46*, 1120-1134.

- 9.- Long, C.M.; Suh, H.H.; Koutrakis, P. Characterization of Indoor Particle Sources Using Continuous Mass and Size Monitors; *J. Air & Waste Management Assoc.* 2000, 50, 1236-1250.
- 10.- Klepeis, N.E.; Ott, W.R.; Repace, J.L. The Effect of Cigar Smoking on Indoor Levels and Particles; *J. Expos. Anal. Environ. Epidemiol.* 1999, *9*, 622-635.
- 11.- Ramachandran, G.; Adgate, J.L.; Hill, N.; Sexton, K.; Pratt, G.C.; Bock, D. Comparison of Short-Term Variations (15-Minute Averages) in Outdoor and Indoor PM2.5 Concentrations; *J. Air & Waste Manage. Assoc.* 2000, *50*, 1157-1166.
- 12.- Nelson, P.R.; Conrad, F.W. Sensitivity Correction for Portable Particulate Monitors for Environmental Tobacco Smoke. Presented at the American Association for Aerosol
- 13. Ogden, M.W.; Heavner, D.L.; Foster, T.L.; Maiolo, K.C.; Cash, S.L.;Richardson, J.D.... Nelson, P.R. Personal Monitoring System for Measuring Environmental Tobacco Smoke Exposure; *Environ. Technol.* 1996, *17*, 239-250.
- 14.- Conner, J.M.; Oldaker, G.B., III; Murphy, J.J. Method for Assessing the Contribution of Environmental Tobacco Smoke to Respirable Suspended Particles in Indoor Environments; *Environ. Technol.* 1990, *11*,189-196.
- 15.- Ogden, M.W.; Maiolo, K.C.; Oldaker, G.B., III; Conrad, F.W. Evaluation of Methods for Estimating the Contribution of ETS to Respirable Suspended Particles. In *Indoor Air '90, Proceedings of the 5th International Conference on Indoor Air Quality and Climate*, 2; Toronto, Ontario, Canada, July 29–August 3, 1990; pp 415-420.
- 16.- Jenkins, R.A.; Palausky, A.; Counts, R.W.; Bayne, C.K.; Dindal, A.B; Guerin, M.R. Exposure to Environmental Tobacco Smoke in Sixteen Cities in the United States as Determined by Personal Breathing Zone Sampling; *J. Expos. Anal. Environ. Epidemiol.* 1996, *6*, 473-502.
- 17.-. Maskarinec, M.P.; Jenkins, R.A.; Counts, R.W.; Dindal, A.B. Determination of Exposure to Environmental Tobacco Smoke in Restaurant and Tavern Workers in One US City; *J. Expos. Anal. Environ. Epidemiol.* 2000, *10*, 36-49.
- 18.- Becquemin, R.H.; Bertholon, G.F.; Attoui, M.; Roy, F.; Roy, M.; Dautzenberg, B.; Particle size in the smoke produced by six different types of cigarettes tobacco. *Rev Mal Resp* 2009; 26: e12-e18